


axtair automorpho® plus

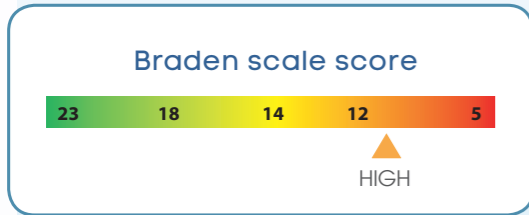

MATTRESS

- ▶ 18 removable and independent Polyurethane Ether cells of 12 cm therapeutic air
 - Head zone: 2 static air cells fitted with a comfort pillow
- ▶ 5 cm polyether foam base located in an independent compartment
- ▶ 4 independent heel relief air cells
- ▶ CPR emergency valve
- ▶ Transport plug: system for opening and sharing 2 air channels (autonomy of 8 hours)
- ▶ Optional cable protection


COVERS AVAILABLE

- ▶ Upper side: stretchable, breathable, impermeable to liquids and permeable to vapour
- ▶ Lower side: non-slip base
- ▶ 2 covers available:
 - PROMUST PU HD (black cover)
 - PROMUST CIC COVER (blue cover) polyamide mesh coated polyurethane/polycarbonate with added silver ions fitted with a non-slip base



PUMP

- ▶ Patented system of automatic and continuous calculation of inflation pressure according to patient's morphology
 - ▶ 3 modes: dynamic, static low pressure, care
 - ▶ Comfort settings: customization of comfort and seating position
 - ▶ Visual and audible alarms.
 - ▶ Indication for servicing
 - ▶ Serial number on the pump's back, simplified instruction and QR code on both sides.

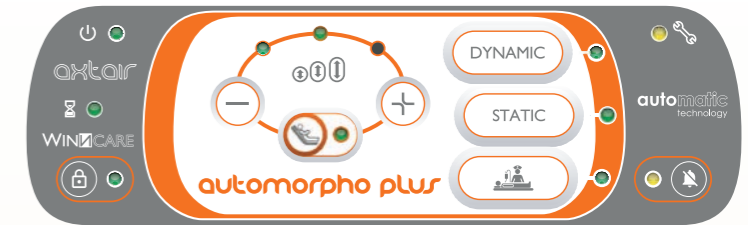

INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.
HELP IN TREATING EXISTING PRESSURE ULCERS FROM STADE 1 TO 4
PATIENT WEIGHT: 30 TO 165 KG/ 4.7 TO 25.9 ST

KEY BENEFITS

- ▶ **Effectiveness and therapeutic safety** for the patient
- ▶ **Easy to use** for caregivers
- ▶ Improved **quality of treatment and patient's quality of life**
- ▶ 3 widths available for **more comfort**
- ▶ **Technical autodiagnostic** by QR code
- ▶ **Ready to use** in a carrying bag
- ▶ **Environmentally-friendly**
- ▶ **Eco-designed and manufactured in France**

DESCRIPTION	REFERENCE	DIMENSIONS / WEIGHT	PATIENT WEIGHT	FIRE STANDARDS	KEY FEATURES
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust PU HD cover For 90 cm wide bed	VAXT4/AUTO-P	195 x 87x17 cm			Overlay mattress CPR emergency valve 18 removable and independent Polyurethane Ether cells of 12 cm therapeutic air 5 cm polyether foam base located in an independent compartment 2 static air cells fitted with a comfort pillow 4 heel reliefs as standard Pump Dimensions: 22 x 25 x 11.5 cm Weight: 2.7 kg Technology based on alternating pressure 1 cell in 2 over a complete cycle from 9 to 14 minutes Patented system of automatic and continuous calculation of inflation pressure Comfort adjustment Visual and audible alarms Flow > 7 l / min Electrical power supply: 220-240 V - 50 Hz Class II type BF Average consumption: 5 Wh Pump decibels: Acoustic pressure according to NF EN3744: <35 dBA Electric and electromagnetic compliance according to the standards IEC 60601-1 Ed.3; IEC 60601-1-2
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust CIC cover For 90 cm wide bed	VAXT4/CIC-P	8kg			
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust PU HD cover For 100 cm wide bed	VAXT4/AUTO100-P	195 x 97x17 cm	30 to 165 Kg	EN 597-1 and 2 GPEM D1 90 Classe D	
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust CIC cover For 100 cm wide bed	VAXT4/CIC100-P	10kg			
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust PU HD cover For 120 cm wide bed	VAXT4/AUTO120-P	195 x 117 x 17 cm			
AXTAIR AUTOMORPHO® PLUS (pump + overlay mattress) Promust CIC cover For 120 cm wide bed	VAXT4/CIC120-P	12kg			

FP000086/08-2018

Cleaning instructions: Promust PU HD and CIC covers



* Behavior of the resistance of the compressor to disinfection with a solution of active chlorine at 1% (10,000 ppm) discontinuous use.

2 years warranty

Photos not contractually binding



WINNOCARE Group - 200 rue Charles Tellier - Actiparc de Grézan - 30034 Nîmes Cedex 1 France
Ph : +33 (0)4 66 02 15 15 - Fax : +33 (0)4 66 02 15 00 - Email : contact.askle@winnocare.fr - www.winnocare.fr

In the interest of constantly improving its products, the WINNOCARE Group reserves the rights to modify the technical specifications of any product without prior notice

Class I medical device in accordance with European Directive 93/42/EEC

